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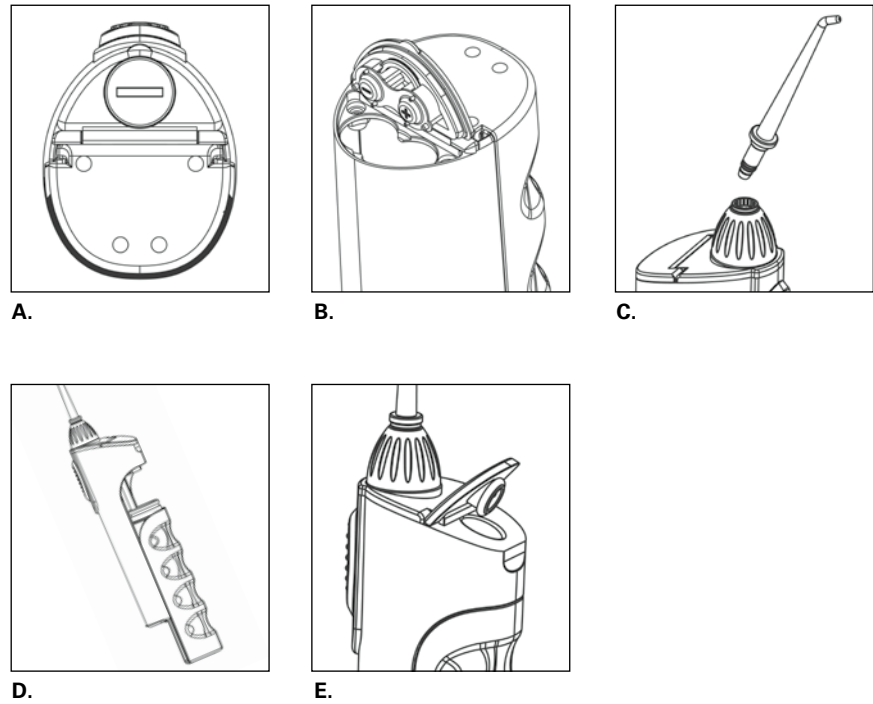
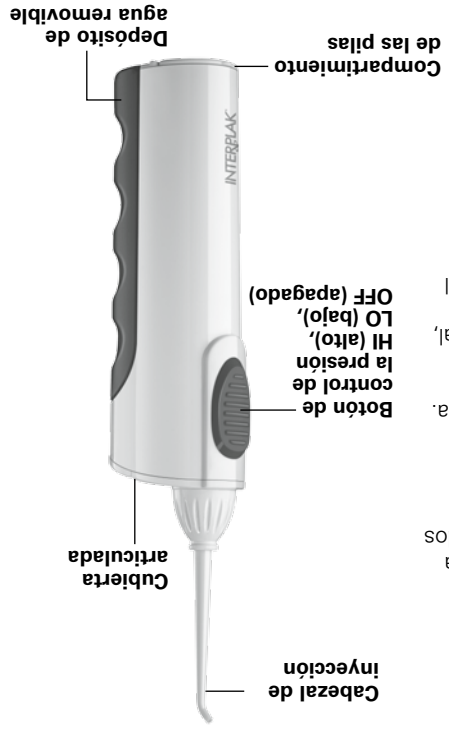
## Instructions for Use:

1. Place water jet tip on unit by snapping tip into place (see C).
  2. Fill water tank. There are two convenient ways to fill the water tank.
    - a. Remove water tank by sliding down the unit (see D). Or,
    - b. Open hinged cover and fill tank while attached to unit (see E).
  - NOTE:** You can also use your favorite dental rinse or mouthwash. Thoroughly rinse reservoir with clean water after each use with mouthwash or dental rinse, and whenever unit has not been used for an extended period of time.
  3. Begin with switch positioned to lowest setting (Lo) to allow gums to become accustomed to water jet cleaning.
  4. Place water jet tip in mouth, and lean over the sink.
  5. Adjust to higher setting (Hi) if desired.
  6. Turn handle to direct the water jet tip toward tooth surfaces and between teeth to rinse away trapped plaque and particles. The pulsing action of the water massages gums as it cleans the entire mouth.
- NOTE:** To prevent damage to delicate gum areas, do not point water jet tip directly at gums. Water jet tip should be held ¼ in. away from your teeth and gums.

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## Initial Preparation:

Before using your Interplak® by Conair® Cordless Portable Water Flossing System for the first time:

1. Turn battery door lock using a coin or flathead screwdriver and insert 2 AA alkaline batteries (see A). Ensure that the batteries are installed correctly with regard to polarity (+ and -) (see B). Replace door. Always purchase the correct size and grade of battery most suitable for intended use. Replace all batteries of a set at the same time.

**NOTE:** Use only alkaline batteries. Do not mix old and new batteries. Do not mix alkaline, standard (carbon-zinc), or rechargeable (nickel cadmium) batteries. Remove batteries from equipment that is not to be used for an extended period of time. When replacing batteries, make sure you dispose of the old batteries properly. Clean the battery contacts and also those of the device prior to battery installation. Remove used batteries promptly.

When replacing the batteries, make sure you dispose of the old batteries in accordance with applicable regulations. They should be disposed of in a location specifically designated for that purpose so they can be recycled safely and without posing any risk for the environment. Do not burn or bury them.


2. Prior to filling the water tank, make sure that the water jet tip is firmly placed onto the top of unit (see C).
3. Allow full reservoir of lukewarm water to pass through the tip before first use. (See Instructions for Use.)

## Welcome to a New Level of At-Home Dental Care

The Interplak® by Conair® Cordless Portable Water Flossing System enhances oral hygiene and significantly contributes to healthy gums when used as part of a responsible oral health management program. It is particularly effective at cleaning around orthodontic appliances, fixed and removable bridge work, implants and splints, and it provides valuable gum massage with every use.

**Remember:** Oral irrigation is a supplement to, not a substitute for, your oral hygiene routine. Ask your dental health care provider to advise you on how the water jet can most effectively meet your specific oral hygiene needs.

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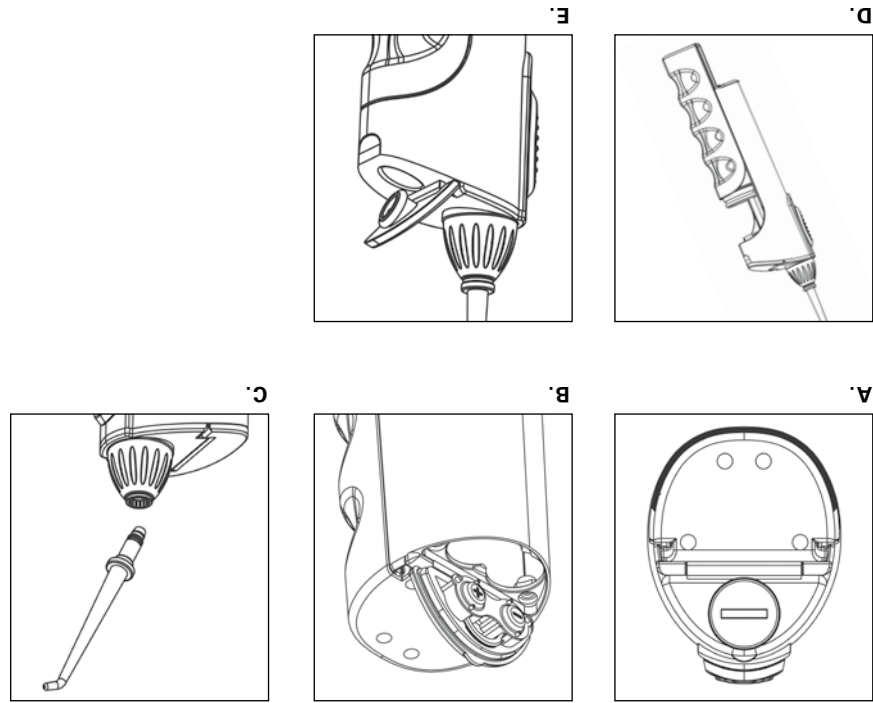**INTERPLAK**  
by CONAIR

*cordless*  
**PORTABLE**  
*water flossing system*

### Instruction Manual • Model WJ4 (All Versions)

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

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### Limited 1-Year Warranty (U.S. and Canada Only)

Conair will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed below, together with your purchase receipt and \$3.00 for postage and handling.

California residents need only provide proof of purchase and should call 1-800-334-4031 for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 12 months from the date of manufacture.

**ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.**

Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

**IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

Please register this product at: [www.conair.com/registration](http://www.conair.com/registration)

### Service Center:

Conair Corporation, Service Department, 7475 N. Glen Harbor Blvd., Glendale, AZ 85307

©2016 Conair Corporation, 150 Milford Road, East Windsor, NJ 08520  
For more information regarding this and other Conair products, visit [www.interplak.com](http://www.interplak.com) or call 1-800-334-4031

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*irrigador bucal inalámbrico*  
**PORTÁTIL**

## IMPORTANT SAFEGUARDS

### READ ALL INSTRUCTIONS BEFORE USING

## DANGER

- DO NOT submerge the unit in water or any other liquids.

## WARNING

- Close supervision is necessary when this product is used by, on, or near children or individuals with certain disabilities.

- Use this product only for its intended use as described in this manual. Do not use attachments unless recommended by the manufacturer.

- Non-rechargeable batteries are not to be recharged. Exhausted batteries are to be removed from the product.

## SAVE THESE INSTRUCTIONS

## SAVE THESE INSTRUCTIONS