



Multi-Speed Wand Massager



Model: UC-569

Operating Manual

Introduction

This Wand Massager is simple to use with 4 speed selections. Powerful vibration massage relaxes tense muscles and relieves pain in your lower back, shoulders, legs or wherever you need it most. Soft rubberized head is suitable for sensitive spots yet firm enough when pressure is needed.

Features of UC-569 Wand Massager:

- Ultra lightweight design
- Soft head, flexible movement for maximum sensation
- Muscle relaxation and pain relief
- 1 Knocking, 3 Vibration speed of up to 6000 rpm

Product Specifications:

Unit size	12.2"×2.4"×2.4"
Unit Weight	1.65lbs
Power Cord Length	98.5"
Alternating Voltage	110-120V
Frequency	60HZ
Vibration level and speed	Level 1: 3500 rpm Level 2: 4200 rpm Level 3: 5000 rpm Level 4: 6000 rpm

Content

Safety Warnings	1~2
Set up	3
Operating the Product	4
Maintenance and Clean.....	5
Trouble shooting	6

Safety Warnings

When the massager is not in use, always unplug from the socket. This prevents electric shock and damage to the product.

Please read these instructions thoroughly and keep for future reference.

Warnings:

Consult your doctor prior to using the product if you have heart, skin and/or other health-related issues.


Please take the following precautions to prevent electric shock, burn, fire and other injuries.

1. Do not leave this product unattended. Unplug from the socket when not in use.
2. Do not use in an excessively humid environment, such as shower, tub, pool or similar.
3. This product is intended for adult use only. Children should be supervised to ensure that they do not play with the appliance.
4. Do not use outdoors.
5. Keep the power cord away from heat sources.
6. Do not disassemble the product unless authorized.
7. This appliance is not intended for use by persons (including children) with reduced physical sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Attention:

1. Do not let babies and children play with the product.
2. Start massaging with the lowest intensity first and gradually go higher if desired.
3. If you feel any discomfort, reduce the intensity or stop using the product immediately.
4. Recommended daily usage is about 20 minutes.
5. Do not use while you are asleep.
6. Do not use on eyes and/or any other sensitive areas
7. To prevent overheating, do not use under blanket, pillow, quilt or similar.
8. Do not use if skin is inflamed, swollen or any other abnormality.

Set up

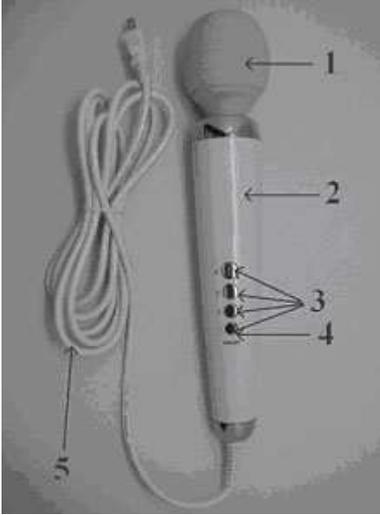

1. Massage head
2. Body
3. Speed selection buttons
4. On/Off button

Operating the product

1. Press the ON/OFF button to turn on the massager, which is at the lowest speed with knocking massage (Level 1)
2. Press Buttons 2-4 to adjust the vibration speed of the massager (Level 4 being the strongest).
3. Massage the body using the sides of the vibrating massage head.
4. Press the ON/OFF button to turn off the massage or unit will automatically turn off after 20 minutes of continuous usage.
5. Recommended time of usage: 20 minutes

Maintenance and Clean

Unplug the product when it needs to be cleaned. Only use alcohol or volatile spray cleaner to clean the product.

1. Make sure that the massager is dried completely prior to use
2. Store massager in a cool and dry place, away from direct sunlight.

Notice: Do not wrap the power cord around the product.

If you follow the operating instructions in this manual, this product will last for a long time.

However, if the following situations happen, please stop using immediately, unplug from socket and contact Sunpentown:

1. If the massage head becomes loose.
2. If the power cord is worn off or the product does not work properly.



Your Guarantee

If this product is found to be faulty as a result of faulty materials or workmanship within one year from date of purchase, it will be repaired free of charge.

This guarantee is subject to the following terms:

- Sunpentown must be notified of the fault.
- Proof of purchase must be presented to Sunpentown's nominated representative.
- The warranty will be void if the product is modified, misused or repaired by an unauthorized person.
- The warranty after repair will not be extended beyond the original one-year period.
- All replacement parts will be new or reconditioned.
- Parts, which are replaced, become the property of Sunpentown.
- The warranty applies for the use of the product in the USA only.

What is NOT COVERED:

- Warranty does not include freight charges.
- Incidental or consequential damage caused by possible defects with this product.
- Damage to product caused by improper power supply voltage, accident, fire, floods or acts of nature.
- Failure of product resulting from unauthorized modifications to the product.
- Improper installation or failure to perform the necessary maintenance.

This GUARANTEE is in addition to your Statutory Rights

SUNPENTOWN INTERNATIONAL INC.
14625 Clark Ave. City of Industry, CA 91745
Tel: 800-330-0388
service@sunpentown.com
www.sunpentown.com